

ottobock.

The new Avantgarde. Living ahead!

Quality for life



At first sight and then forever!

Of course it is love, since love is where it all begins. Only those who love themselves will also be loved, since self-worth has a lot to do with feelings. And you are worth it. For you, we have once again significantly improved the best way to get around in an active wheelchair. For you, we have combined design and function in the perfect symbiosis of aesthetics and technology. For you, we have built the new Avantgarde. We hope you will fall in love with the new Avantgarde. Ideally at first sight.

The new Avantgarde is available in custom colour and material combinations as well as three versions: The CV designed for maximum variability, the CS for sporty dynamics and the ultra-lightweight CLT. Show your true colours with the new Avantgarde.

Love for details Avantgarde



Being different means thinking differently

Those of us who are not dependent on a wheelchair to get around have to learn how to think like a wheelchair user. This is far from simple. Which is why we got help from people who know what they are talking about for the development of the new, third generation of the Avantgarde. We examined and questioned every detail, every function, every material, every screw and every seam in the course of intensive development efforts. We tested, discarded, improved and tested again. Until there was nothing left that we could improve.

We learned a lot from the wheelchair-bound members of the Otto Bock development team. Especially how important it is to not merely settle for being different, but to learn how to think differently in order to develop a new, self-confident self. Which is precisely why we built the new Avantgarde.



Avantgarde CV in Snake Clear Red, equipped with anti-tipper, crutch holder, elevating leg rest right, swing-away side panels with long armrests, quick-release axles with tetra-unlock, height-adjustable push handles and wheel lock lever extension.

Avantgarde CV

The all-round talent

Perfect handling

Every day serves up fresh challenges. With the new Avantgarde CV, you can look forward to them – because you receive precisely the support you desire with sporty performance, ergonomic handling and outstanding functionality. The special details include a vertical anterior frame and swing-away, removable footrest holders allowing you to approach furniture, the bed or bathtub with precision while also facilitating a straightforward transfer. The footrest holders can be unlocked even with limited hand functionality.

Basic technical data

Load capacity	140 kg
Weight	approx. 12 kg

Transport weight

Frame	approx. 7.800 g
Footrest	approx. 900 g
Drive wheel	from 1.250 g



Avantgarde CS

Winning image

Sporty design in Pure Metal Green

The new Avantgarde CS is ideal for all those who want to ride through life boldly, letting nothing slow them down. Thanks to outstanding agility and manoeuvrability, the folding frame active wheelchair is equally convincing in confined spaces and over long distances. Special features of the "Powerbox" include compact exterior dimensions with a highly stable design and closed anterior frame with integrated footrest. Frame angles of 70 and 80 degrees allow the legs to be positioned at the desired angle. A 90-degree frame without footrest for bilateral amputees rounds out the selection.

Basic technical data	
Load capacity	140 kg
Weight	approx. 11 kg
Transport weight	
Frame	approx. 8.500 g
Drive wheel	from 1.250 g





The new Avantgarde CS in Pure Metal Green, offset frame in 80°, equipped with Froglegs Softrolls and Froglegs suspension forks, shock absorber system, SPOX rear wheels, knee lever wheel locks with short operating lever, tube armrests, business seat upholstery and "Backster" fixed back shell.



The new Avantgarde CLT in Pure Metal Blue, permanently welded, frame in 70°, equipped with 4" Froglegs Softrolls, Infinity Ultralight drive wheels, carbon clothing protector side panels, high-frequency welded back upholstery made of 3D spacer fabric and welded push handles.

Avantgarde CLT

The lightweight

9,200 grams of ergonomics

Demonstrating individuality has never been easier: At just 9,200 grams, the CLT is the lightweight of the new Avantgarde series. This significant weight advantage of the folding frame active wheelchair is achieved through the modular lightweight design and tubular aluminium frame with central double crossbrace. A drive wheel adapter supports fine-tuning of the tipping point, wheelbase, seat angle, wheel camber and toe-in. For experienced users who are confident in their sitting position, a permanently welded version without adapters is available as an alternative offering an even further weight reduction.

Basic technical data	
Load capacity	100 kg
Weight adjustable	approx. 10 kg
Weight welded	approx. 9 kg
Transport weight	
Frame adjustable	approx. 8.000 g
Frame welded	approx. 7.000 g
Drive wheel	from 1.250 g



Frame colours

Customising your Avantgarde

Design the new Avantgarde according to your wishes with the colour and modular add-on components of your choice. Precisely meeting your day-to-day requirements is possible.

The new Avantgarde by Otto Bock Mobility Solutions: Wheelchairs – made exactly for you.

Pure metal red




Candy red



Pure metal bronze



Sparkle light orange

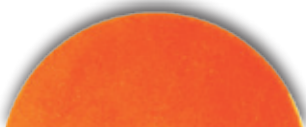

Snake clear red



Dormant vinho sparkle



Sparkle granny smith



Pure metal green



Dormant marine blue



Dormant sky blue



Pure metal blue



Ice blue




White pearl



Brilliant silver




Anthracite metallic

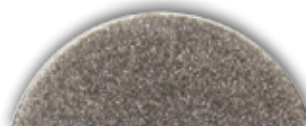

Deep black matte



Clear reptile glossy

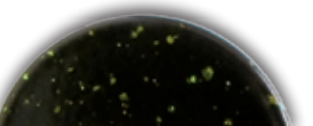

Listen to your heart. And follow your mind. Talk to your Ottobock partner at your medical supply company, health centre or specialised rehabilitation dealer and find out together which Avantgarde is the best match for you. We are sure you will find it. We know it will be exactly right.

Further information is now also available under www.avantgarde3.com

My Avantgarde Not later, but now





Otto Bock HealthCare GmbH
Max-Näder-Straße 15, 37115 Duderstadt
T +49 5527 848-1590 · F +49 5527 848-1676
reha-export@ottobock.de · www.ottobock.com